

01-18-00

A

Please type a plus sign (+) inside this box →

+

PTO/SB/05 (12/97) (modified)  
Approved for use through 09/30/00 OMB 0651-0032  
Patent and Trademark Office U.S. DEPARTMENT OF COMMERCE

01/14/00

UTILITY  
PATENT APPLICATION  
TRANSMITTAL

(Only for new nonprovisional applications under 37 CFR 1.53(b))

Attorney Docket No. LABV.P0002

First Named Inventor: Tseung

Title: Method and Apparatus of Automatic Tissue Staining

Express Mail Label No. EL497530733US

## CERTIFICATE OF MAILING BY "EXPRESS MAIL"

Express Mail Label No. EL497530733US

Date of Deposit: November 3, 1999

I hereby certify that this paper or fee is being deposited with the United States Postal Service "Express Mail Post Office to Addressee" service under 37 C.F.R. § 1.10 on the date indicated above and is addressed to: Assistant Commissioner for Patents, Washington, D.C. 20231

Dag Johansen

## APPLICATION ELEMENTS

See MPEP chapter 600 concerning utility patent application contents

## ADDRESS TO:

Assistant Commissioner for Patents  
Box Patent Application  
Washington, DC 20231

1. ☐ Fee Transmittal Form  
(Submit an original, and a duplicate for fee processing)
2. ☒ Specification -- Total Pages 27 --  
Including  
- Title page  
- Background of the Invention  
- Brief Summary of the Invention  
- Brief Description of the Drawings  
- Detailed Description  
- Claim(s)  
- Abstract of the Disclosure
3. ☒ Drawing(s) (35 USC 113) -- Total Sheets 18 --
4. ☐ Oath or Declaration (unsigned) -- Total Pages 1 --  
a. ☐ Newly executed (original or copy)  
b. ☐ Copy from a prior application (37 CFR 1.63(d))  
(for continuation/divisional with Box 17 completed)  
[Note Box 5 below]  
c. ☐ DELETION OF INVENTOR(S)  
Signed statement attached deleting inventor(s) named in  
the prior application, see 37 CFR 1.63(d)(2) and 1.33(b)
5. ☐ Incorporation By Reference (useable if Box 4b is checked)  
The entire disclosure of the prior application, from which a copy of the  
oath or declaration is supplied under Box 4b, is considered as being  
part of the disclosure of the accompanying application and is hereby  
incorporated by reference therein

6. ☐ Microfiche Computer Program (Appendix)
7. ☒ Nucleotide and/or Amino Acid Sequence Submission  
(if applicable, all necessary)  
a. ☐ Computer Readable Copy  
b. ☐ Paper Copy (identical to computer copy)  
c. ☐ Statement verifying identity of above copies

## ACCOMPANYING APPLICATION PARTS


8. ☐ Assignment Papers (cover sheet & document(s))
9. ☐ 37 CFR 3.73(b) Statement ☐ Power of Attorney  
(when there is an assignee)
10. ☐ English Translation Document (if applicable)
11. ☐ Information Disclosure Statement (IDS)/PTO-1449 ☐ Copies of IDS Citations
12. ☐ Preliminary Amendment
13. ☒ Return Receipt Postcard (MPEP 503)  
(Should be specifically itemized)
14. ☐ Small Entity ☐ Statement filed in prior application,  
Statement(s) Status still proper and desired
15. ☐ Certified Copy of Priority Document(s) (if foreign priority is claimed)
16. ☐

17. If a CONTINUING APPLICATION, check appropriate box and supply the requisite information:

☐ Continuation ☐ Divisional ☐ Continuation-in-part (CIP) of prior application No. \*

## 18. CORRESPONDENCE ADDRESS

Dag Johansen, Esq.  
STATTLER JOHANSEN & ADELI LLP  
P.O. Box 51860  
Palo Alto, California 94303-0728  
Telephone: (650) 934.0470 ext. 102  
Facsimile: (650) 934-0475

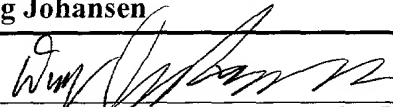Name	Dag Johansen	Registration Number:	36,172
Signature		Date:	January 14, 2000

APPLICATION FOR UNITED STATES LETTERS PATENT

FOR

**Method and Apparatus for Automatic Tissue Staining**

Inventors: KEN K. TSEUNG  
GLENN TAKAYAMA  
NORMAN K. RHETT  
MARK V. CORL

Prepared by:  
Stattler, Johansen & Adeli  
P.O. Box 51860  
Palo Alto, California 94303-0728  
Phone: 650-934-0470  
Fax: 650-934-0475

"Express Mail" mailing label number EL497530733US

Date of Deposit 1/14/2000

I hereby certify that this paper or fee is being deposited with the United States Postal Service "Express Mail Post Office to Addressee" service under 37 CFR 1.10 on the date indicated above and is addressed to the Commissioner of Patents and Trademarks, Washington, D.C. 20231.

Dag Johansen  
(Typed or printed name of person mailing paper or fee)

Dag Johansen 1/14/2000  
(Signature of person mailing paper or fee)

## Method and Apparatus for Automatic Tissue Staining

### FIELD OF THE INVENTION

5                   The present invention relates to the field of medical lab equipment. In particular the present invention discloses a fully automated system for staining tissue specimens and cell preparations.

### BACKGROUND OF THE INVENTION

10                   A normal function of medical laboratories is to examine cells and cell tissue under a microscope. The lack of contrast between individual cells and the background matrix or between individual parts of cells can make it difficult to examine cell and tissue preparations. To improve the contrast, researchers have applied stains to  
15 cell and tissue specimens to be examined. The stains are absorbed differently by the various structures in cells such that the contrast between the different cell structures is improved.

                  Staining tissue specimens is a nontrivial time consuming process. Often a  
20 number of different staining and rinsing stages are required. Each stage requires a specific amount of reagent or buffer and takes a specific amount of time. Thus, trained technicians are often employed to perform such operations. Furthermore, hospitals and laboratories must stain very large numbers of tissue specimens for patient diagnoses. Thus, automated tissue staining systems have been developed. By automating the  
25 process, expensive human labor is eliminated. Furthermore, automatically staining

specimens significantly reduces the probability of an error occurring during the staining process.

- To ensure that the proper staining procedures are followed, most
- 5 automatic staining systems require that the user carefully enter the staining protocol and load the proper reagents. The complicated procedures require user training before such devices can be operated effectively. It would therefore be desirable to simplify the operation of an automatic tissue-staining device.

## SUMMARY OF THE INVENTION

A method and apparatus for specimen slide preparation is disclosed. The method and apparatus of the present invention uses slide trays that have receptacles for at least one specimen slide and an associated reagent pack. The specimen slide and/or reagent pack includes an identifier that specifies a particular slide preparation protocol that should be followed. The method and apparatus reads the identifier to determine the particular slide preparation protocol and then prepares the specimen slide according to the particular slide preparation protocol. The apparatus may obtain some or all of the reagents needed for the particular slide preparation protocol from the reagent pack.

Other objects feature and advantages of present invention will be apparent from the drawings and the detailed description that follows.

## BRIEF DESCRIPTION OF THE DRAWINGS

The objects, features, and advantages of the present invention will be apparent to one skilled in the art, in view of the following detailed description in which:

5

**Figure 1A** illustrates a cut-away view of a first tiltable sink assembly embodiment.

**Figure 1B** illustrates the tiltable sink assembly of **Figure 1A** wherein the tiltable sink assembly is tilted to the right.

10

**Figure 1C** illustrates the tiltable sink assembly of **Figure 1A** wherein the tiltable sink assembly is tilted to the left.

15

**Figure 2A** illustrates a cut-away view of a second tiltable sink assembly embodiment.

**Figure 2B** illustrates the tiltable sink assembly of **Figure 2A** wherein the tiltable sink assembly is tilted to the right.

20

**Figure 2C** illustrates the tiltable sink assembly of **Figure 2A** wherein the tiltable sink assembly is tilted to the left.

**Figure 3** illustrates a conceptual diagram of the internal fluid flow components of the autostainer apparatus.

25

**Figure 4a** illustrates a first arrangement of a reagent pack for the autostainer apparatus of the present invention.

5                   **Figure 4b** illustrates a second arrangement of a reagent pack for the autostainer apparatus of the present invention.

**Figure 5a** illustrates a four well reagent pack or the autostainer apparatus of the present invention.

10                   **Figure 5a** illustrates a six well reagent pack or the autostainer apparatus of the present invention.

15                   **Figure 5c** illustrates an eight well reagent pack or the autostainer apparatus of the present invention.

**Figure 6a** illustrates a first embodiment of bulk packaged reagent packs.

20                   **Figure 6b** illustrates a second embodiment of bulk packaged reagent packs.

**Figure 7a** illustrates a front view of a combined slide and reagent rack for preparing four slides.

**Figure 7b** illustrates a back view of a combined slide and reagent rack for preparing four slides.

**Figure 8A** illustrates a flow diagram of the general operation of the  
5 autostainer apparatus.

**Figure 8B** illustrates a flow diagram summary of the slide protocol scheduling system of the autostainer control program.

10 **Figure 9a** through **9d** illustrate a flow diagram describing the details of adding a new slide to the Currently Scheduled Array (CSA).



## DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENT

A method and apparatus for automatically staining tissue specimens is disclosed. In the following description, for purposes of explanation, specific  
5 nomenclature is set forth to provide a thorough understanding of the present invention. However, it will be apparent to one skilled in the art that these specific details are not required in order to practice the present invention. For example, the present invention has been described with reference to staining of tissue specimens. However, the same techniques can easily be applied to other types of slide preparation work.

### **The Autostainer Hardware**

The present invention comprises advances in automated slide staining. An example of an automatic slide staining apparatus can be found the first figure of U.S.  
15 Patent 5,839,091, issued November 17, 1998, and entitled "Method and apparatus for automatic tissue staining" which is hereby incorporated by reference in its entirety. An autostainer system is used for staining tissue specimens that are placed onto glass slides. The present invention uses several different slide racks wherein each slide rack holds one or more slides. In one embodiment, there are six slide racks and each slide rack is  
20 capable of holding four slides such that the autostainer can prepare twenty-four different slides simultaneously.

An autostainer uses a robotic delivery system that delivers bulk reagents, small supply reagents, buffer solutions, and air to the glass slides. The robotic delivery  
25 system is controlled by a computer system. The computer system executes an autostainer

control program that sends control commands to control the robotic delivery system. In one embodiment, the robotic delivery system of the autostainer consists of an X-axis mechanism, a Y-axis mechanism, and a Z head as illustrated in U.S. Patent 5,839,091. The Z head has a Bulk fluid dispensing tube for dispensing a few selected bulk reagents and buffer rinse solution, an air blade to blow air onto slides, and a syringe probe for picking up reagents that will be placed onto the glass slides.

To prevent contamination, the syringe probe is cleaned in a reagent probe wash bin between the uses of different reagents. The wash bin has three different receptacles that are used in three stages. The first hole is used to rinse the inside of the probe by forcing buffer rinse solution through the inside of probe and down into a first drain receptacle. The second receptacle is used to clean the outside of the probe by forcing buffer rinse solution through the inside of probe while the probe is in the tightly confined second receptacle such that the buffer solution is force upward on the outside of the probe. Finally, the probe is placed into a third receptacle and air is forced through the probe to clean out the buffer rinse solution.

Beneath the slide racks of the autostainer is a sink assembly. The sink assembly catches the reagents and buffer rinse solution that drip off the slides. **Figure 1A** illustrates a cut-away top view of a first embodiment of a sink assembly. As illustrated in **Figure 1A**, a tiltable sink **210** sits beneath a support bracket **205** for the slide racks (not shown). The tiltable sink **210** may be tilted left or right using a tilt mechanism **220**. When the tiltable sink **210** is tilted down on the right side as illustrate in **Figure 1B**, all the liquid waste spilling off the slides will drain out of drain hole **231** on the right side, through corrugated tubing **241** and drain pipe **251**, and finally to a waste

line #1 (not shown). Similarly, when the tiltable sink **210** is tilted down on the left side as illustrated in **Figure 1C**, all the liquid waste will drain out of drain hole **232**, through corrugated tubing **242** and drain pipe **252**, and finally to a waste line #2 (not shown). As illustrated in **Figures 1A** through **1C**, the two waste lines conveniently exit the device at the same location.

With the tiltable sink system of the present invention, the present invention can send different types of down different waste lines such that a first waste line may be used to remove nonhazardous waste and the other waste line may be used to remove hazardous waste. The non hazardous waste line may simply be connected to a sewage drain pipe. The hazardous waste line should be connected to hazardous waste container that is disposed of appropriately.

**Figures 2A** through **2C** illustrate an alternate embodiment of a tiltable sink system. Referring to **Figures 2A** through **2C**, the drain line travels straight down such that the two waste lines exit at different locations.

#### Fluid Flow components

Several other components are also located inside the autostainer. **Figure 3** illustrates a conceptual diagram the fluid flow components of the autostainer. Referring to **Figure 3**, three different output devices located on the Z-Head assembly of the autostainer deliver air or fluid to the slides.

A first output device is the bulk fluid dispensing tube **321**. The bulk fluid dispensing tube **321** is used to deliver relatively large quantities of fluid to slides. The

bulk fluid dispensing tube **321** may dispense buffer solution from buffer supply **301** or reagents from internal bulk reagent supplies **302**, **303**, and **304**. The particular fluid dispensed is selected by 8-way distribution valve **331**.

5                   A second output device is the reagent-dispensing probe **322**. The reagent-dispensing probe **322** draws in specific reagents using syringe pump **340** and then dispenses drawn-in reagents onto specific slides. To prevent contamination of the reagent-dispensing probe **322** from different reagents, the autostainer cleans the reagent-dispensing probe **322**. Specifically, the autostainer first flushes the reagent dispensing probe **322** by forcing buffer solution through the reagent dispensing probe **322** using  
10                   three-way valve **332**. Then, the autostainer dries the reagent dispensing probe **322** by forcing air through the reagent dispensing probe **322** using three-way valve **333**.

                    The third output device is air blade **323**. Air blade **323** is used to dry off  
15                   slides and blow away extra reagents.

#### Autostainer Reagent Packs

                    To simplify operation, the present invention introduces an autostainer system with greatly simplified operation. To use the autostainer system of the present  
20                   invention, a user simply adds a set of slides where each slide is accompanied by a specific reagent pack that contains the reagents needed for a specific slide preparation protocol. The reagent pack further includes information that identifies the slide preparation protocol to be performed.

**Figure 4a** illustrates a first embodiment of a reagent pack that may be used with the autostainer system of the present invention. As illustrated in **Figure 4a**, the reagent pack comprises a container with a set of wells (**401** to **406** in **Figure 4a**) for storing reagents. The container is sealed using a cover **410**. The cover **410** may include an identification mark such as barcode **411** to identify the slide preparation protocol to be performed. The cover **410** protects and retains the reagents in the wells **401** to **406**, yet the reagent dispensing probe **322** is able to puncture and access the reagents as needed.

The cover **410** may further include a second barcode identification sticker **420**. The barcode identification sticker **420** can be placed directly onto a slide to be performed such that the autostainer will automatically know the slide preparation protocol to be performed.

**Figure 4b** illustrates an alternate embodiment of a reagent pack **490**. In the alternate embodiment of **Figure 4b**, the identifiers **431** and **441** have been placed parallel to the line of wells **401** to **406**. In yet another embodiment (not shown), a tiltable sink is created with four drain holes with one hole at each corner. In such an embodiment, four different waste systems may be used.

The autostainer system of the present invention is capable of handling many different slide preparation protocols that require different numbers of reagents. **Figures 5a, 5b, and 5c** illustrate different sized reagent packs for different protocols. **Figure 5a** illustrates a simply four reagent pack for simple slide preparations. **Figure 5b** illustrates the six-reagent pack of **Figure 4a**. **Figure 5c** illustrates an eight-reagent pack for complex slide preparation protocols.

The autostainer will be used to perform large numbers of slide preparations. Each slide preparation requires a reagent pack. To simplify purchasing, reagent packs may be purchased in bulk packages. **Figure 6a** illustrates one possible bulk package where several reagent packs are sold in a perforated two-dimensional matrix of individual reagent packs. **Figure 6b** illustrates an alternate bulk package where several reagent packs are sold as a strip of connected of individual reagent packs.

#### Autostainer Slide Racks

The slide racks for the autostainer of the present invention have been designed for ease-of-use and maximum flexibility. **Figure 7a** illustrates a front view of one possible slide rack that may be used in the autostainer of the present invention. The slide rack **700** has four slide positions such that the slide rack **700** can hold four slides **710**. Other slide rack embodiments may hold more or less than four slides. The slides may be standard U.S. or international sized slides. The slide rack **700** further includes four reagent positions for storing reagent packs **720** associated with the adjacent slide positions. **Figure 7b** illustrates a back view of the slide rack **700**.

#### **Autostainer Control and Programming**

As stated in the previous section, the autostainer is controlled by a computer system. In a present embodiment, the computer system is based on a standard Personal Computer (PC) motherboard with a PCI bus. The computer system runs an autostainer control program to control the operation of the autostainer.

The autostainer control program is a sophisticated control program that implements many security, automatic slide protocol programming, control, and logging features. To fully describe the autostainer control program, this document will step through a sample use of the autostainer.

5

#### User Loading

To operate the autostainer, a user loads slide trays (as illustrated in **Figures 7a** and **7b**) with slide specimens and associated reagent packs. In one embodiment, the user simply places the proper reagent pack in the reagent pack receptacle adjacent to the specimen slide to be prepared. In another embodiment, the user places an associated identifier on the frosted area of the specimen slide to be prepared. This can be done by peeling off a barcode sticker **420** from the reagent pack **400** and placing it on the frosted area of the slide. By specifically labeling each slide with the protocol to be performed, no error can be introduced by accidentally placing the wrong reagent pack next to a specimen slide.

10

15

#### Slide Preparation Protocol Identification

After loading one or more slide trays, the user places the loaded slide trays into the autostainer system. The user then instructs the autostainer system to commence slide preparation by activating a “restart” input. If the user has added very high priority slides to the autostainer, the user may instead press a “STAT” input to indicate that the new slides are very high priority. The autostainer system then commences operation by first examining the loaded slides and reagent packs to determine the slide preparation protocols that need to be performed. Specifically, the autostainer system reads all the identifiers (on the slides and/or reagent packs) and then consults a slide preparation

20

25

protocol database that maps the different identifiers on the reagent packs with the slide preparation protocols to be performed.

- 5 The slide preparation protocol database may be periodically updated by shipping media such as floppy disks or CD-ROMs that are inserted into an appropriate drive on the autostainer. In one embodiment, the autostainer system uses a network to automatically retrieve database updates. This may be performed by having the autostainer coupled to a telephone line with a modem such that the autostainer automatically calls specific number to obtain slide preparation protocol database updates.
- 10 In an alternate embodiment, the autostainer may be coupled to the global Internet such that the autostainer may connect to a server that stores the most recent slide preparation protocol database information. In any of these embodiments, the autostainer may also simultaneously receive program updates such that the software that controls the autostainer may be automatically updated.

15

#### Slide Staining Apparatus Operation

- Once a user has loaded the autostaining apparatus, the user may start staining run. **Figure 8A** illustrates the general procedure of operating the autostaining device. Referring to **Figure 8A**, the user loads a slide/stain tray at step **801**. The user
- 20 then starts the apparatus by pressing the restart or stat button at step **802**. The autostaining apparatus then examines all the slides to determine the slide protocols that must be performed at step **803**. The system may compare a barcode sticker **402** on a slide with a barcode **411** on the cover **410** of an adjacent reagent pack to ensure that the proper reagent pack has been placed next to each slide. After examining all the slides and



reagent packs, the system creates a staining schedule. Details on how the staining apparatus creates the staining schedule will be presented in the following section.

5 After calculating a staining schedule, the autostaining apparatus beings to process slides according to a created staining schedule. The system continues the staining operations until one of the conditions of step **805** are detected (or an error occurs). Specifically, step **805** determines if the staining run is completely done, the staining run is partially done, or if the user has pressed the pause button.

10 If the staining run is partially done, then the system proceeds to step **807** where the user is informed of the completed slide(s). The user may wish to remove those slides such that they can be examined.

15 If the user pressed the pause button, then the system temporarily ceases operation. The user may then load additional slide trays at step **801**. The system remains paused as it waits for the user to press the restart or stat button. Once the user has pressed either the restart or stat button at step **802**, the system proceeds to step **803** to create a new staining schedule. The system then resumes slide staining at step **804** using the new staining schedule.

20 Referring back to step **805**, if the system is completely done with the staining run, then the system may cease operation. As shown in **Figure 8A**, a user may continually remove completed slide trays and add new slide trays continually by only using the pause and restart/stat buttons. Thus, the staining apparatus may be in  
25 continuous operation all day long.

### Slide Preparation Scheduling

Once the slide preparation protocols that need to be performed have been identified, the autostainer control program proceeds to calculate the most efficient  
5 dispensing pattern for performing the desired slide protocols.

**Figure 8B** illustrates a summary of how the autostainer control program calculates the most efficient dispensing pattern. Referring to **Figure 8B**, at step **810** the autostainer control program first determines if a run was already in progress. If a run was  
10 already in progress, then the autostainer control program adds the current slides into a scheduling table at step **815**. The scheduling table will be used to program the autostainer control program.

Next, at step **820**, the autostainer control program has the z-head assembly  
15 scan all the slide positions to determine if any new slides have been added or any existing slides have been taken away. If a slide that had not yet completed its protocol was removed, the autostainer may inform the user of the error and ask if the user wishes to place the slide back into the autostainer. When one or more new slides have been added, the autostainer control program adds those slides into the current scheduling table as  
20 stated in step **820** of **Figure 8B**. If the user designated the new slides as “STAT”, then those newly added slides are given top priority.

After the autostainer control program has added the filled the scheduling table, the autostainer control program proceeds to step **830** to begin automatic  
25 programming. At step **830**, the autostainer control program first places the slides into a

specific order. In one embodiment, the autostainer control program orders the slides using this order:

1. Highest priority slides (such as STAT slides).
2. Slides with the longest incubation time.
3. Numerical order.

The system then proceeds to re-order the slides using a well-defined method. **Figure 8B** illustrates one possible method of ordering the slides. Initially, the first slide from the order of step **830** defines the currently scheduled array (CSA). The CSA defines a newly proposed slide order. At step **845**, the method selects the next unscheduled slide from the scheduling table. That slide is then placed into the CSA as stated in step **850**. At step **860**, the method determines if the slide fit properly into the current stage of the CSA. The method of determining if a slide “fits” into the current stage of the CSA is fully defined in **Figures 9a to 9d**. If the slide did not fit, the method proceeds to step **870** to determine if the control program has tried to fit all the remaining unscheduled slides into the current stage. If there are unscheduled slides that the autostainer control program has not attempted to put into the current stage, then the method proceeds to step **873** to select an unscheduled slide with a different protocol. The autostainer control program then returns to step **850** to attempt to add that slide into the current stage.

If, at step **870**, the autostainer control program is unable to fit any unscheduled slide into the current stage, then the autostainer control program move to the next stage of the current schedule array (CSA) as stated in step **879**. The autostainer control program then proceeds back to step **845** to start placing slides into the new stage.

Referring back to step **860**, after each new slide is added to the CSA, the autostainer control program determines if all the slides have been scheduled as stated in

step **880**. If all the slides have been scheduled, the autostainer control program determines if it has created a slide-processing schedule that is faster than the current best time. If the autostainer control program determines that it created a slide-processing schedule faster than the previous best time, then the autostainer control program saves the newly created schedule as the best schedule at step **883**. At step **890**, the autostainer control program determines if it has tried all the possible slide orders. If it has not, the autostainer control program proceeds to step **895** where the autostainer control program changes the slide order while maintaining the designated slide priority order. The autostainer control program then proceeds back to step **840** to test a new slide schedule.

5 If, at step **890**, the autostainer control program determines that it has tried every possible slide order (that maintains the slide priority order), the autostainer control program returns the best schedule.

10

The foregoing has described a method and apparatus for automatic tissue and cell preparation staining. It is contemplated that changes and modifications may be made by one of ordinary skill in the art, to the materials and arrangements of elements of the present invention without departing from the scope of the invention.

15

## CLAIMS

We claim:

- 1                    1.        A method of operating an autostainer device, said method  
2 comprising the steps of:  
3                    accepting a slide tray, said slide tray having at least one specimen slide and a  
4                    reagent pack associated with said specimen slide, said reagent pack  
5                    comprising a first identifier that specifies a particular slide preparation  
6                    protocol;  
7                    reading said first identifier from said reagent pack; and  
8                    preparing said specimen slide according to said particular slide preparation  
9                    protocol.
- 1                    2.        The method as claimed in claim 1 wherein said reagent pack is  
2 associated with said specimen slide by being adjacent to said specimen slide.
- 1                    3.        The method as claimed in claim 1 wherein said reagent pack is  
2 associated with said specimen slide by having a second identifier on said specimen slide  
3 that is the same as said first identifier.

1                   4.       The method as claimed in claim 1 wherein said reagent pack  
2 comprises a set of wells, each well containing a reagent needed for said particular slide  
3 preparation protocol.

1                   5.       The method as claimed in claim 1 wherein said reagent pack  
2 comprises a peel-off identifier, said peel-off identifier for placement on said specimen  
3 slide.

1                   6.       A method of operating an autostainer device, said method  
2 comprising the steps of:  
3                   accepting a slide tray, said slide tray having at least one specimen slide and a  
4                   reagent pack associated with said specimen slide, said specimen slide  
5                   comprising a first identifier that specifies a particular slide preparation  
6                   protocol for said specimen slide;  
7                   reading said first identifier; and  
8                   preparing said specimen slide according to said particular slide preparation  
9                   protocol.

1                   7.       The method as claimed in claim 6 wherein said reagent pack is  
2 associated with said specimen slide by being adjacent to said specimen slide.

1                   8.       The method as claimed in claim 6 wherein said reagent pack is  
2 associated with said specimen slide by having a second identifier that is the same as said  
3 first identifier.

1                   9.       The method as claimed in claim 6 wherein said reagent pack  
2 comprises a set of wells, each well containing a reagent needed for said particular slide  
3 preparation protocol.

1                   10.      The method as claimed in claim 6 wherein said reagent pack  
2 comprises a peel-off identifier containing said first identifier, said peel-off identifier for  
3 placement on said specimen slide.

1                   11.      An apparatus for staining specimen slides, said apparatus  
2 comprising:  
3           more than one slide tray, said slide tray for holding more than one specimen slide;  
4           an automatic staining head assembly, said automatic staining head assembly for  
5           depositing reagents on said specimen slides, said automatic staining head  
6           assembly further comprising an input device for reading identifiers that  
7           specify slide preparation protocols to perform;

8 a control system, said control system coupled to said automatic staining head  
9 assembly for controlling said automatic staining head assembly to prepare said  
10 specimen slides during a staining run;  
11 a pause input, said pause input for pausing said apparatus during said staining run;  
12 and  
13 a restart input, said restart input for restarting said apparatus after adding new  
14 specimen slides onto on of said slide trays;  
15 wherein said control system causes said automatic staining head assembly to read a new  
16 set of identifiers associated with said new specimen slides to add said new specimen  
17 slides to said staining run.

1 12. The apparatus as claimed in claim 11 wherein said apparatus  
2 further comprises reagent packs.

1 13. The apparatus as claimed in claim 12 wherein each said reagent  
2 pack comprises a set of wells, each well containing a reagent needed for said particular  
3 slide preparation protocol.

1 14. The apparatus as claimed in claim 12 wherein said identifiers  
2 comprise a barcode on each said reagent pack.



1                   15.    The apparatus as claimed in claim 11 wherein said identifiers  
2   comprise a set of barcodes on said specimen slides.

1                   16.    The apparatus as claimed in claim 11 further comprising:  
2   a STAT restart input, said STAT restart input for restarting said apparatus after  
3   adding new specimen slides onto on of said slide trays wherein said new  
4   specimen slides are given high priority;

1                   17.    A slide rack for a slide staining apparatus, said slide rack  
2   comprising:  
3   a first receptacle for accepting a specimen slide; and  
4   a second receptacle for accepting a reagent pack, said reagent pack containing at  
5   least one reagent needed to prepare said specimen slide.

1                   18.    The slide rack as claimed in claim 17 further wherein said reagent  
2   pack further comprises an identifier that identifies a slide preparation protocol for said  
3   specimen slide.

1                   19.    The slide rack as claimed in claim 17 further wherein said first  
2   receptacle and said second receptacle are adjacent to each other.

1                    20.     A reagent pack for a slide staining apparatus, said reagent pack  
2 comprising:  
3                    a set of wells, said well containing reagents for a specific slide preparation  
4                    protocol; and  
5                    an identifier, said identifier associated with said slide preparation protocol.

1                    21.     The reagent pack as claimed in claim 20 wherein said identifier  
2 comprises a peel-off sticker for placement on an associated specimen slide.

1                    22.     The reagent pack as claimed in claim 20 further wherein said  
2 identifier comprises a peel-off sticker for placement on an associated specimen slide.

1                    23.     A slide staining apparatus, said apparatus comprising:  
2 at least one slide rack for holding a slide specimen to be prepared; and  
3 a tiltable sink assembly, said tiltable sink assembly having a first drain hole on a  
4 first side such that liquid material drains through said first drain hole when  
5 tilted down on said first side, said tiltable sink assembly having a second drain  
6 hole on a second side such that liquid material drains through said second  
7 drain hole when tilted down on said second side.

1                    24.    The apparatus as claimed in claim 23 wherein said first drain hole  
2    is coupled to a sewage system.

1                    25.    The apparatus as claimed in claim 23 wherein said second drain  
2    hole is coupled to a hazardous waste container.


1                    26.    The apparatus as claimed in claim 23 wherein said second drain  
2    hole is coupled to a corrugated tube.

1                    27.    An apparatus for staining specimen slides, said apparatus  
2    comprising:  
3                    at least one slide tray, said slide tray for holding at least one specimen slide and  
4                    an associated reagent pack, said associated reagent pack having reagents  
5                    needed for processing said specimen slide; and  
6                    an automatic staining head assembly, said automatic staining head assembly for  
7                    obtaining said reagents from said associated reagent pack and depositing  
8                    reagents on said specimen slide.

## Abstract of the Disclosure

A method and apparatus for specimen slide preparation is disclosed. The method and apparatus of the present invention uses slide trays that have receptacles for at least one specimen slide and an associated reagent pack. The specimen slide and/or reagent pack includes an identifier that specifies a particular slide preparation protocol that should be followed. The method and apparatus reads the identifier to determine the particular slide preparation protocol and then prepares the specimen slide according to the particular slide preparation protocol. The apparatus may obtain some or all of the reagents needed for the particular slide preparation protocol from the reagent pack.

# Figure 1A

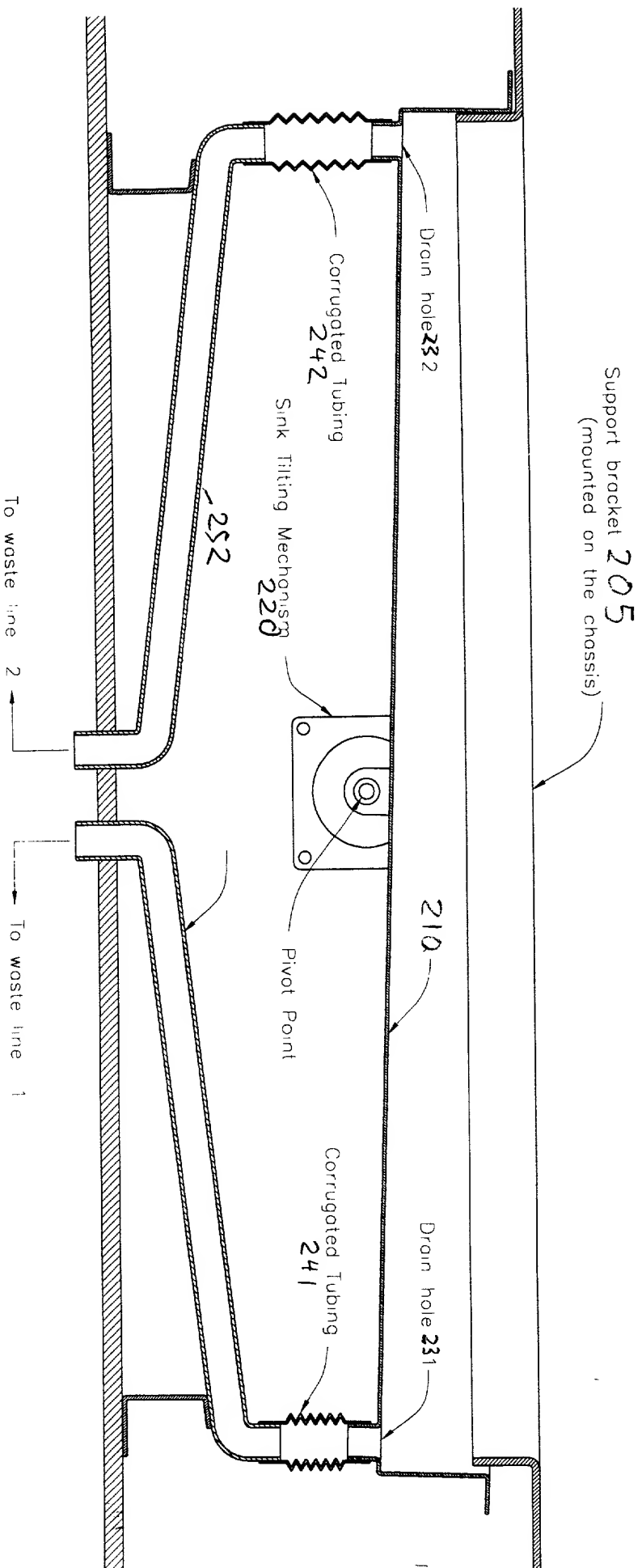

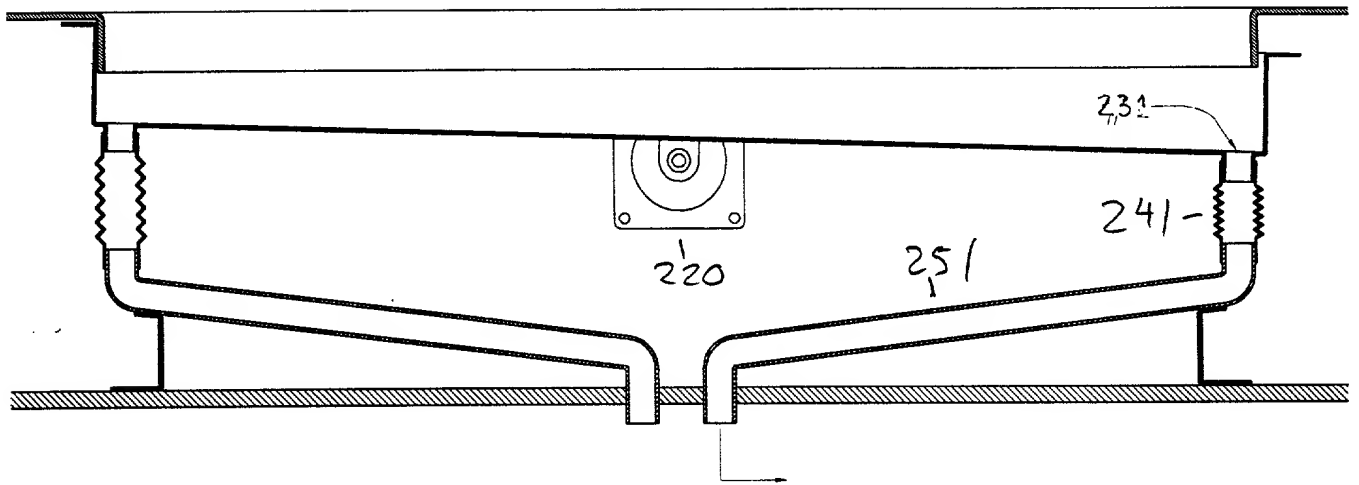


Figure 1B

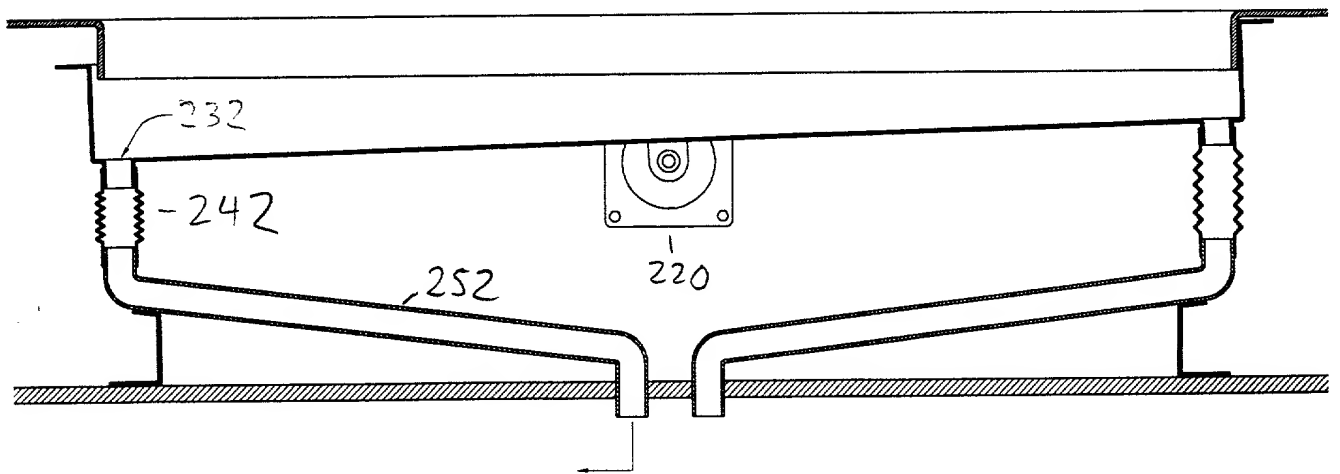
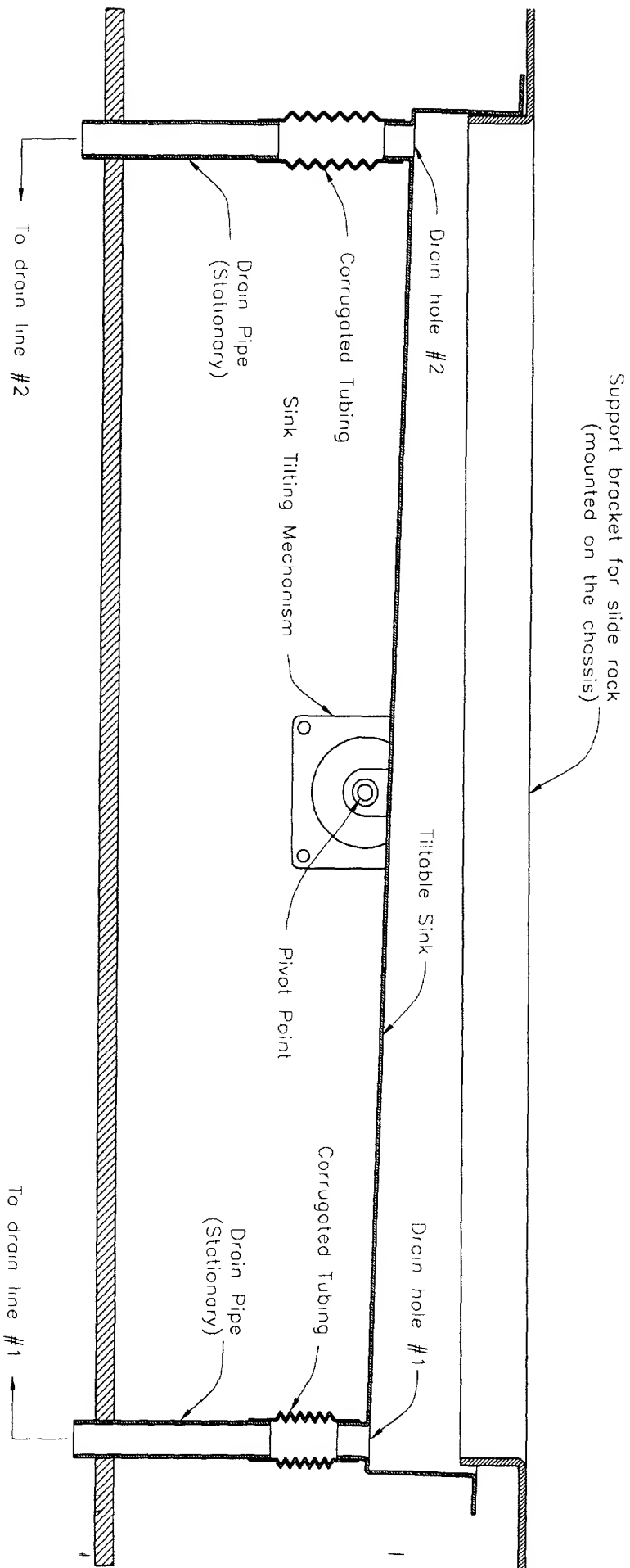


Figure 1C

Figure 2A



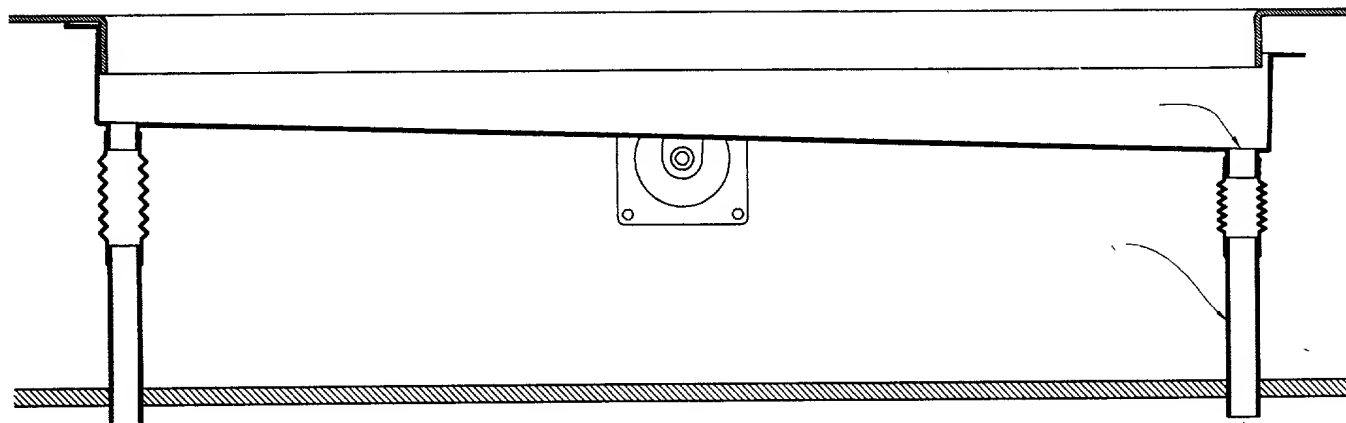


Figure 2B

To drain line #1


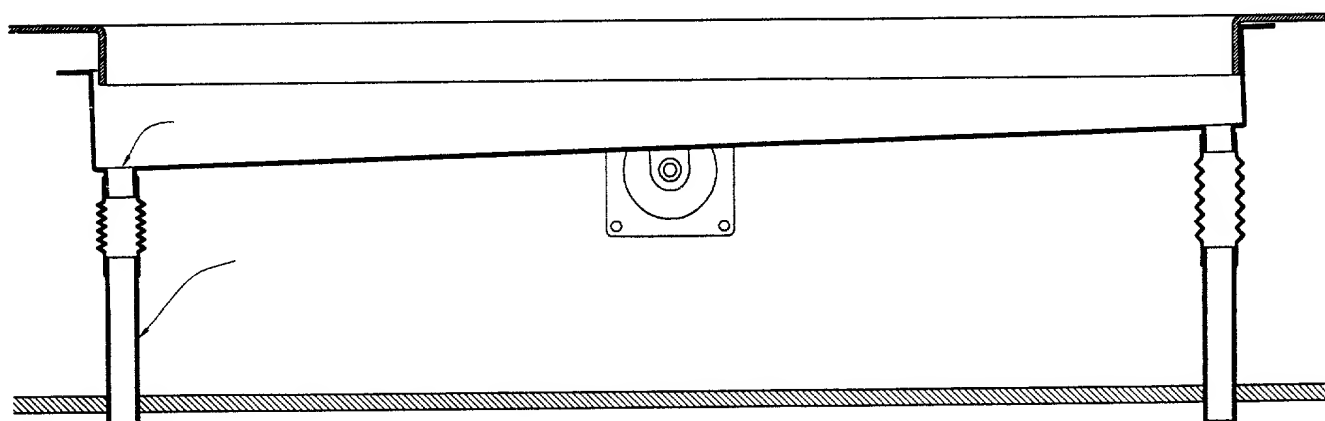


Figure 2C


To drain line #2




# Figure 3

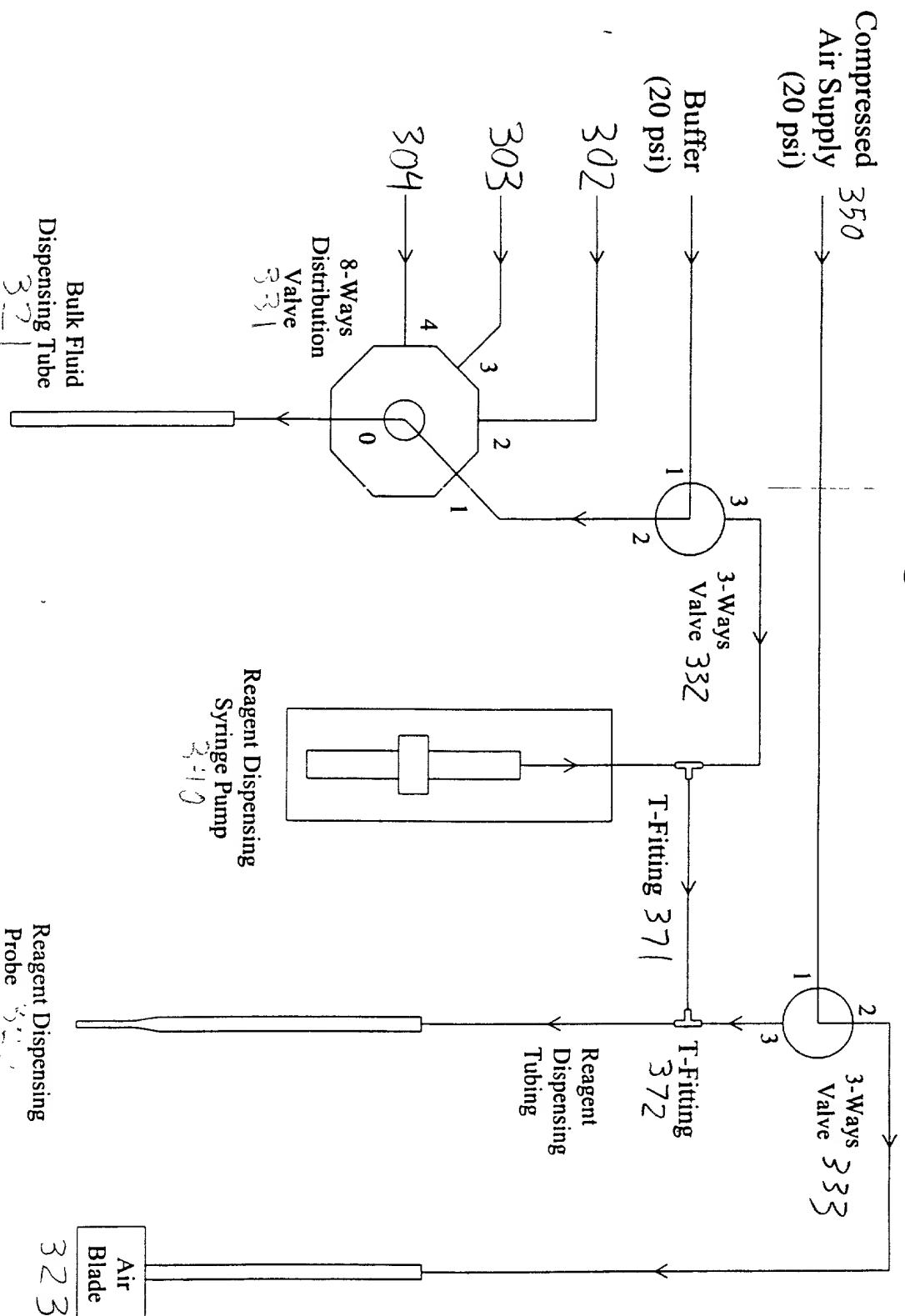


**Figure 4A**

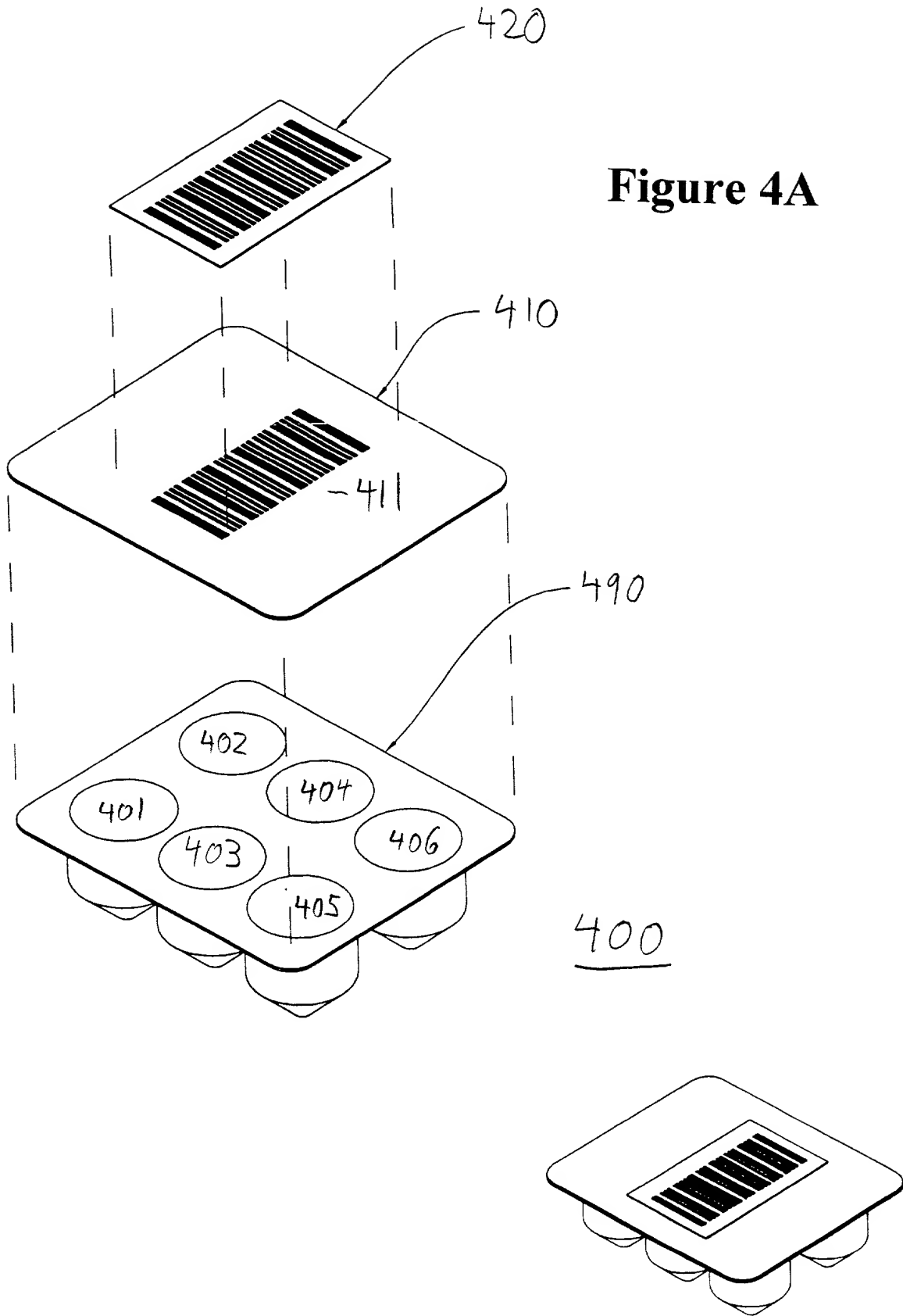

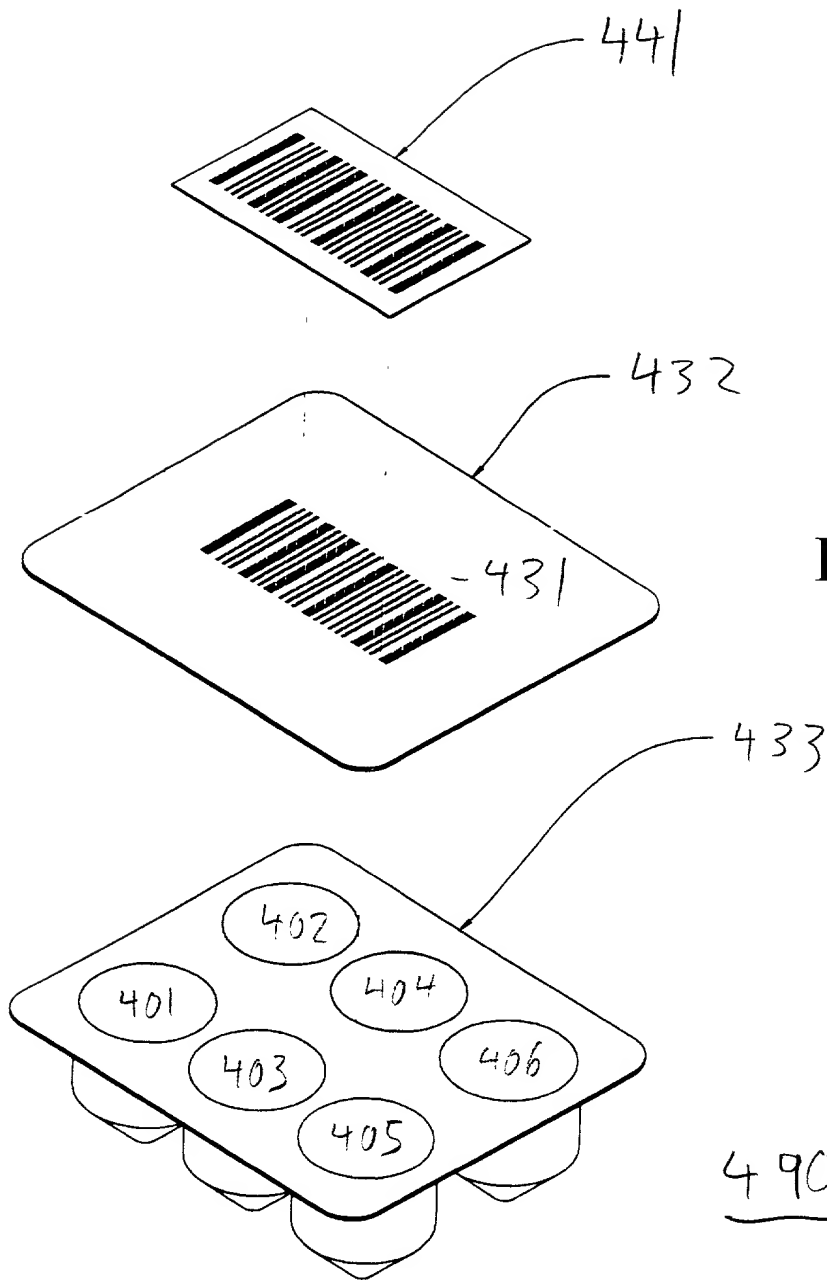


**Figure 4B**

**Figure 5A**

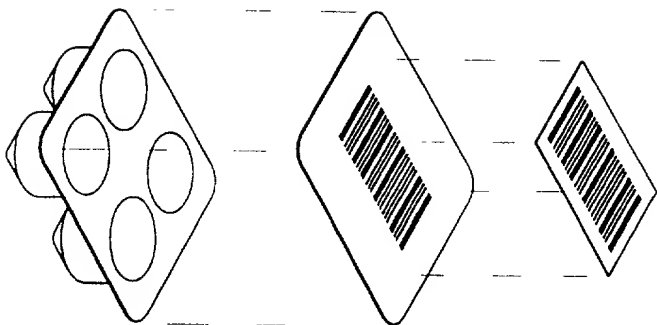


**4 WELLS  
REAGENT PACK**

**Figure 5B**

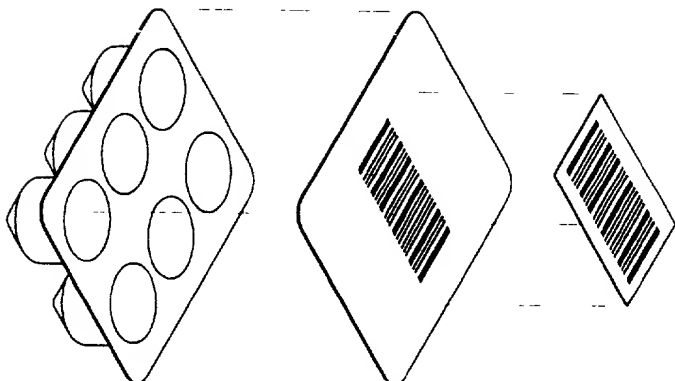


**6 WELLS  
REAGENT PACK**

**Figure 5C**

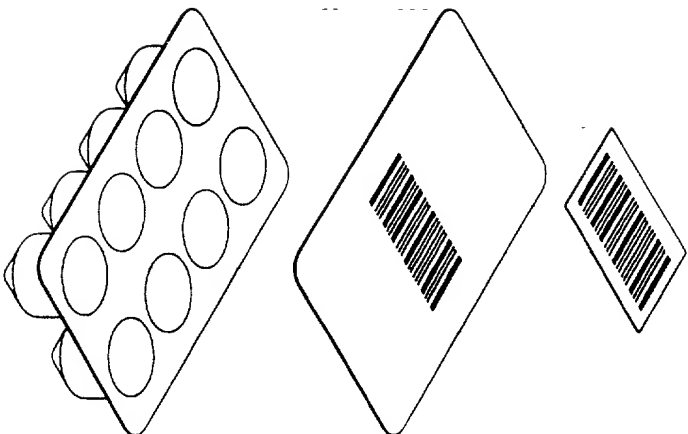

**8 WELLS  
REAGENT PACK**

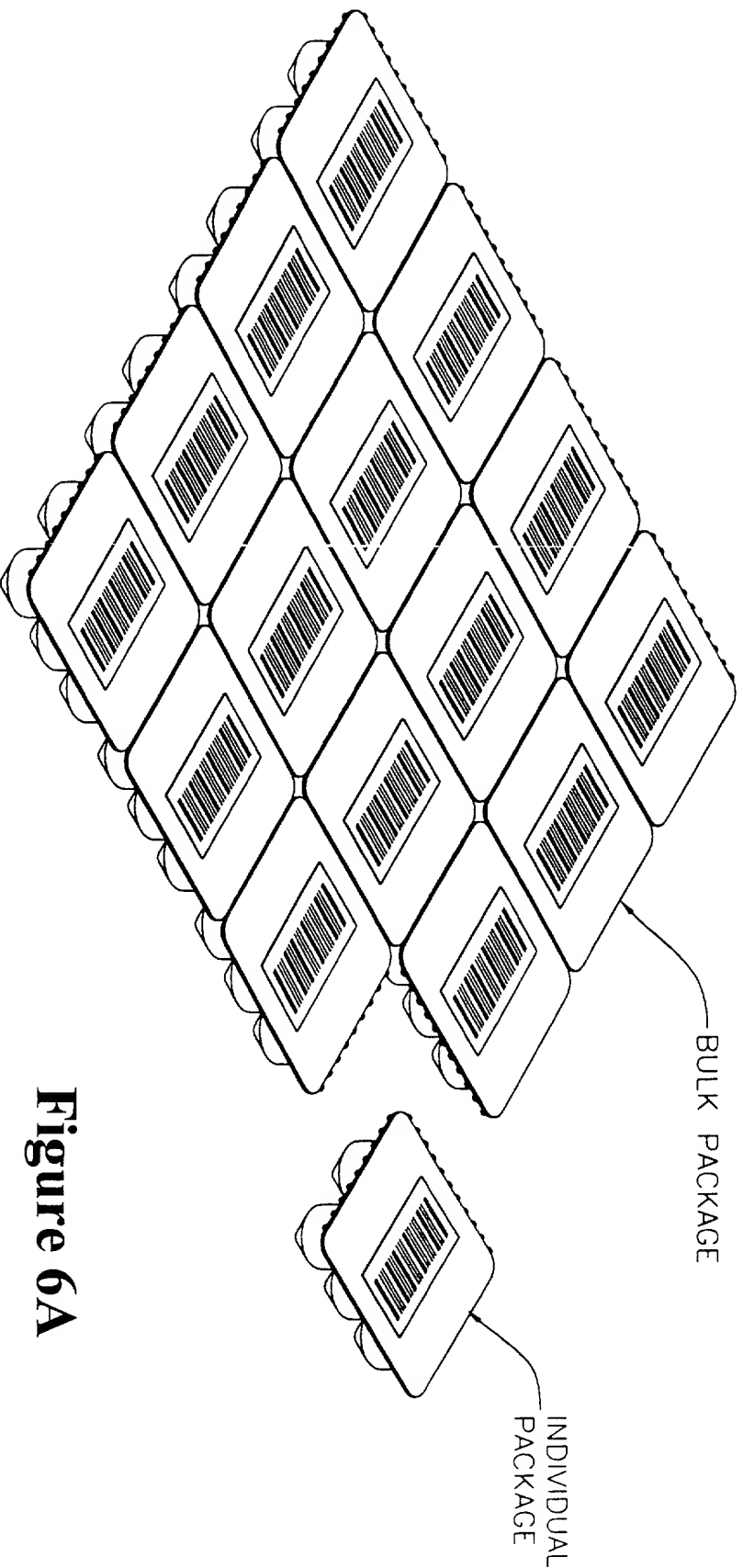

**Figure 6A**


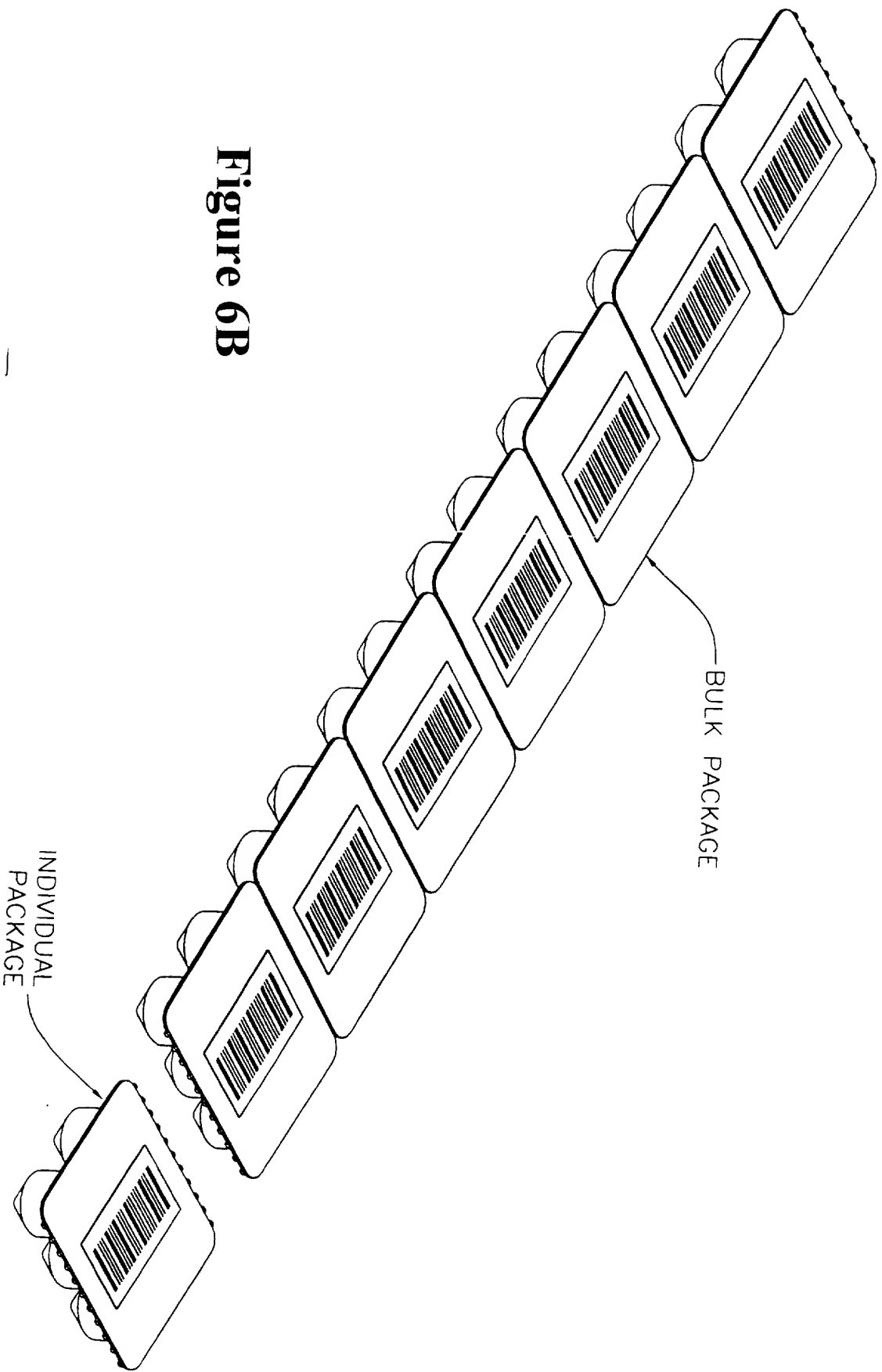



Figure 6B

FIG. 6B is a perspective view of a row of eight rectangular packages, each with a barcode. A label "BULK PACKAGE" points to the middle of the row, and a label "INDIVIDUAL PACKAGE" points to the package at the far right.

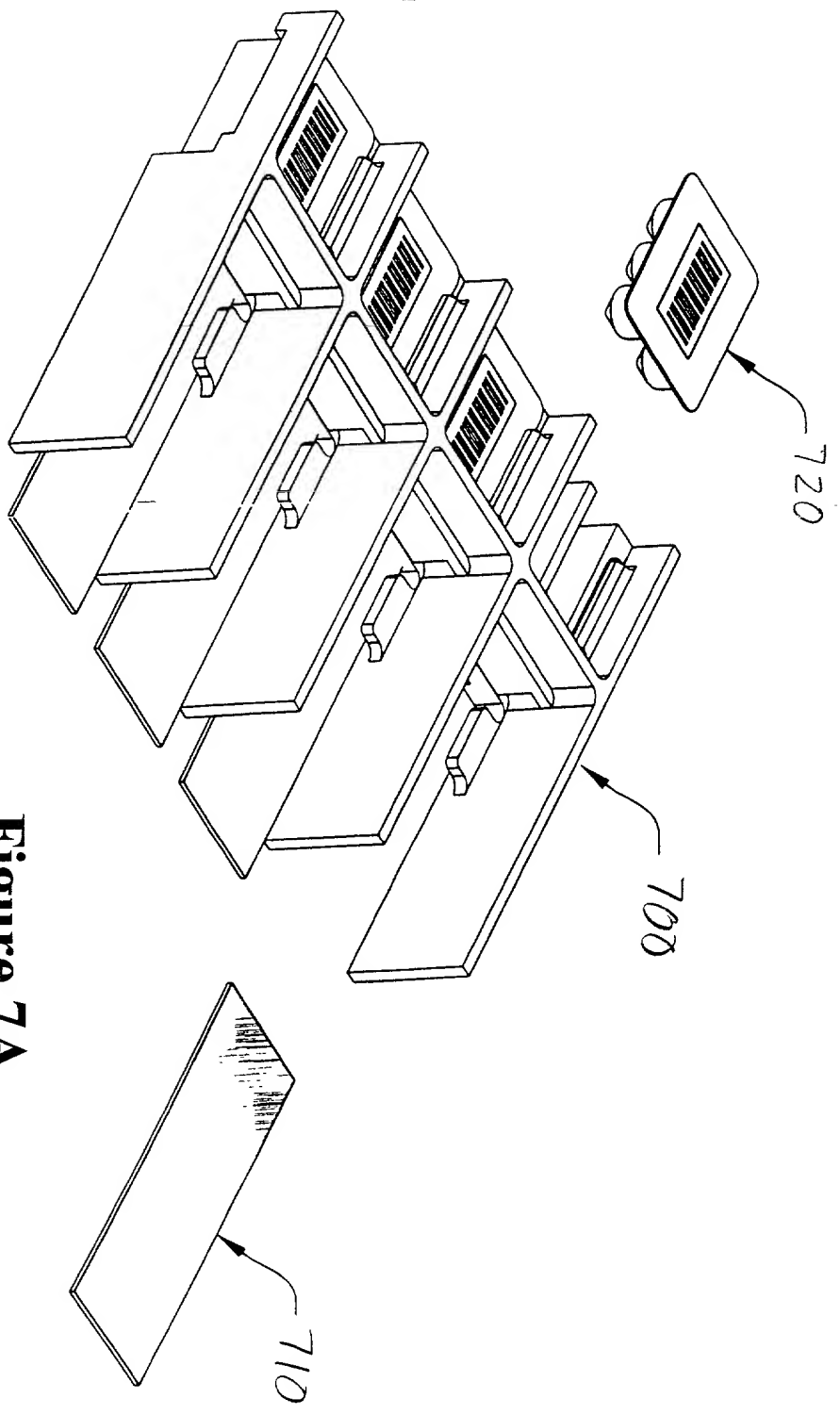

**Figure 7A**


FIG. 7A is a perspective view of a multi-layered assembly 700, 710, 720. The assembly 700, 710, 720 includes a top layer 700, a middle layer 710, and a bottom layer 720. The top layer 700 is a thin, flat plate. The middle layer 710 is a thicker plate with a textured surface. The bottom layer 720 is a small, square component with a grid-like pattern on its top surface.

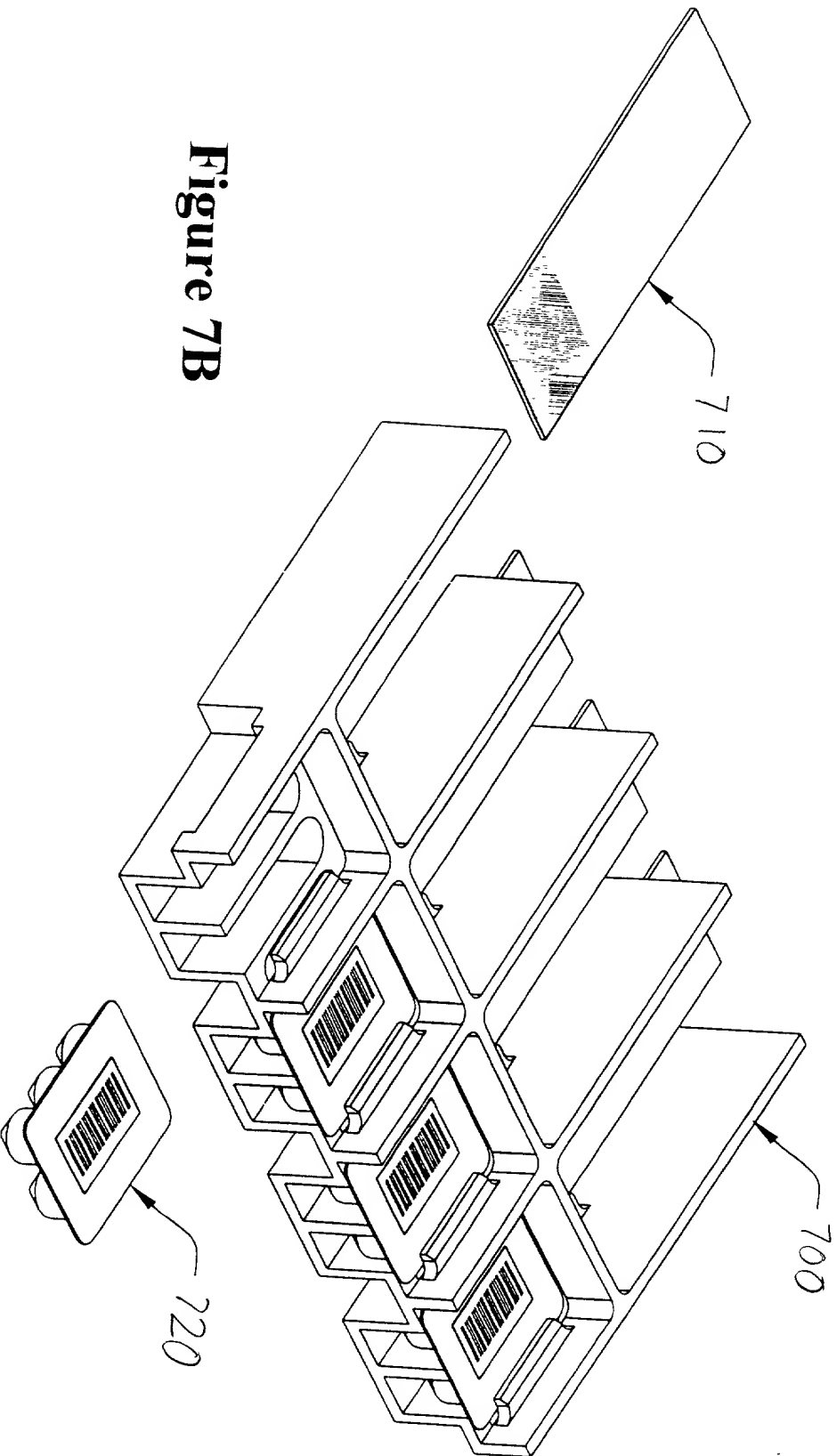


**Figure 7B**



**Figure 8A**

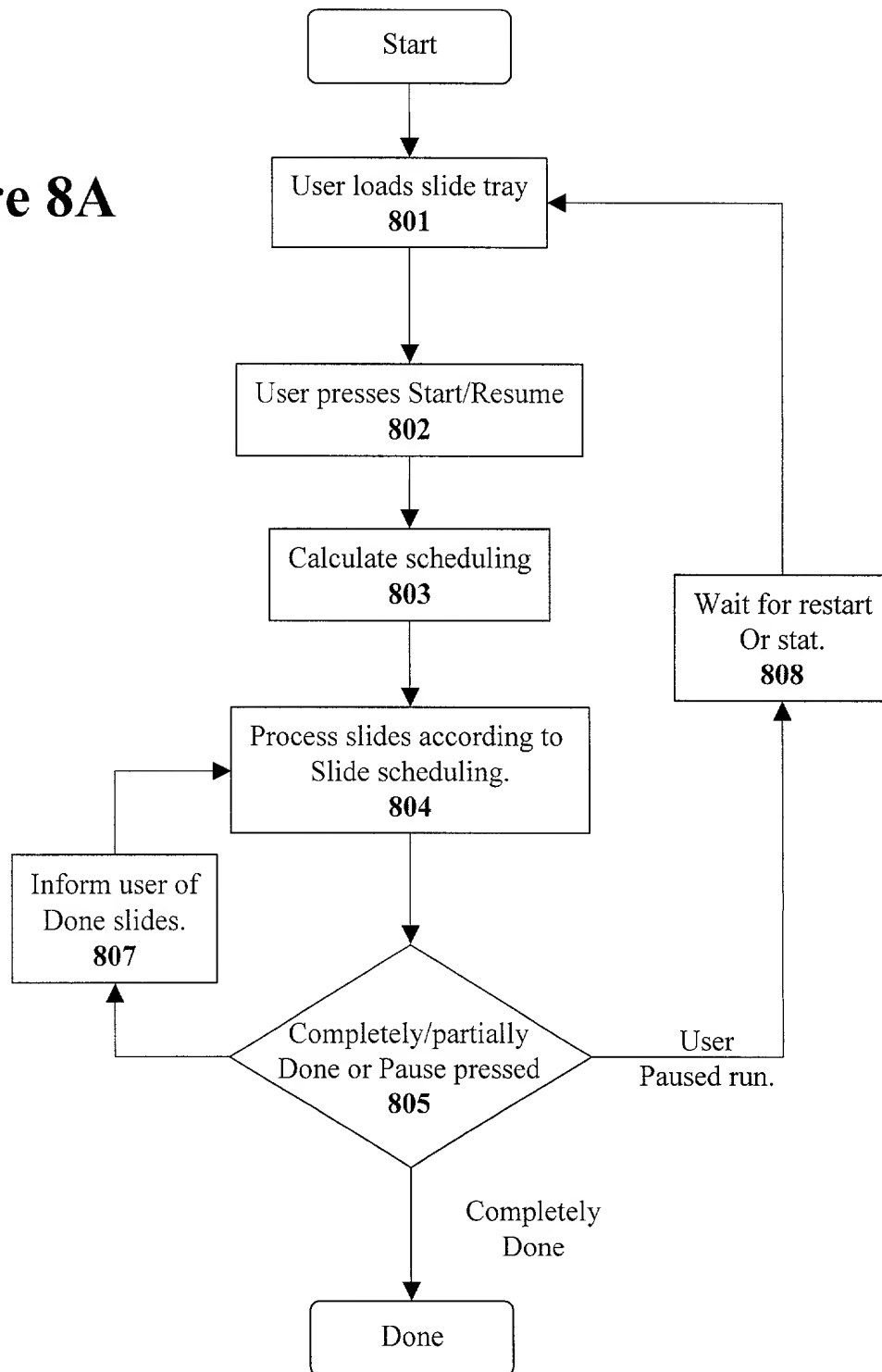

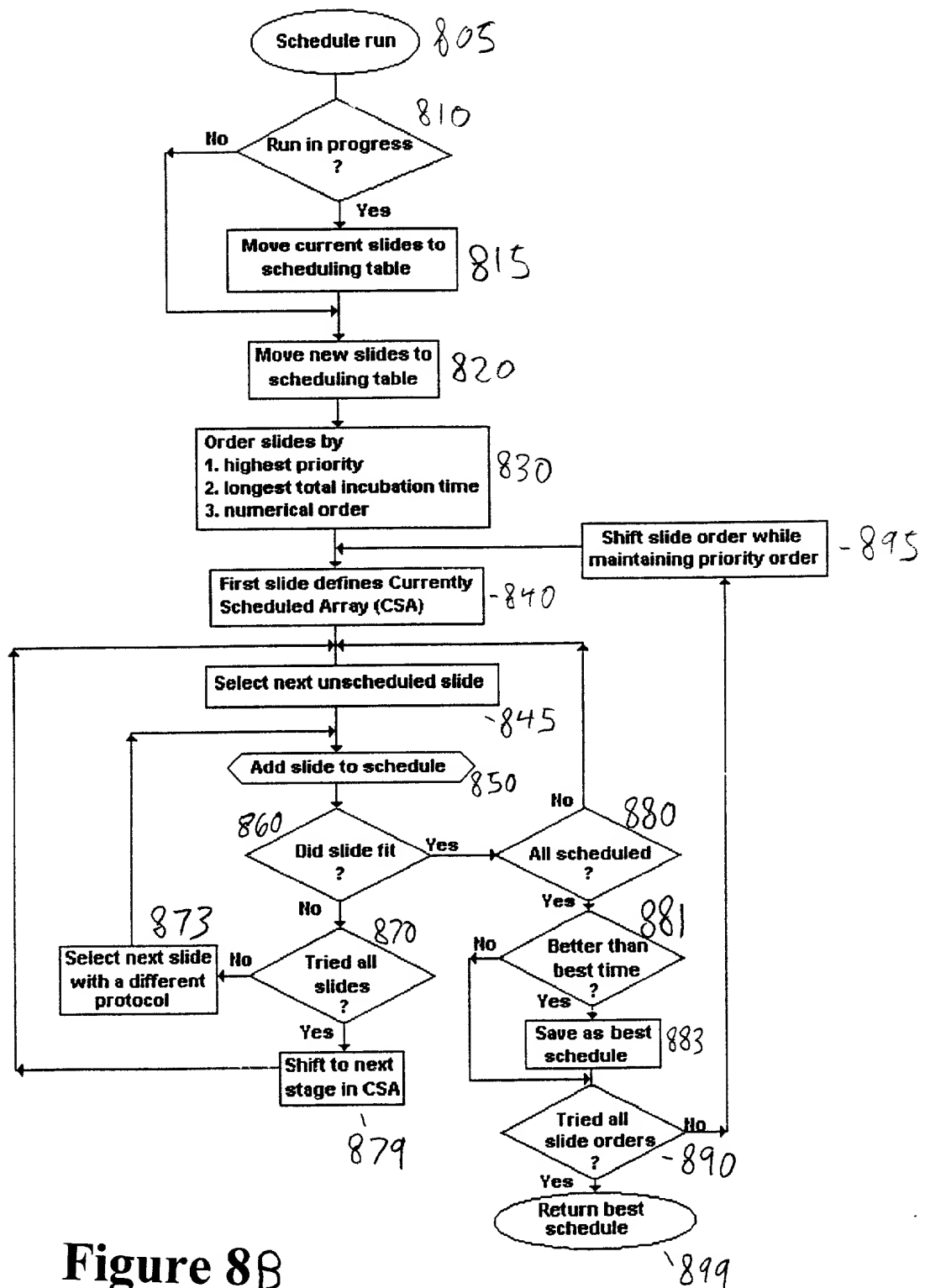

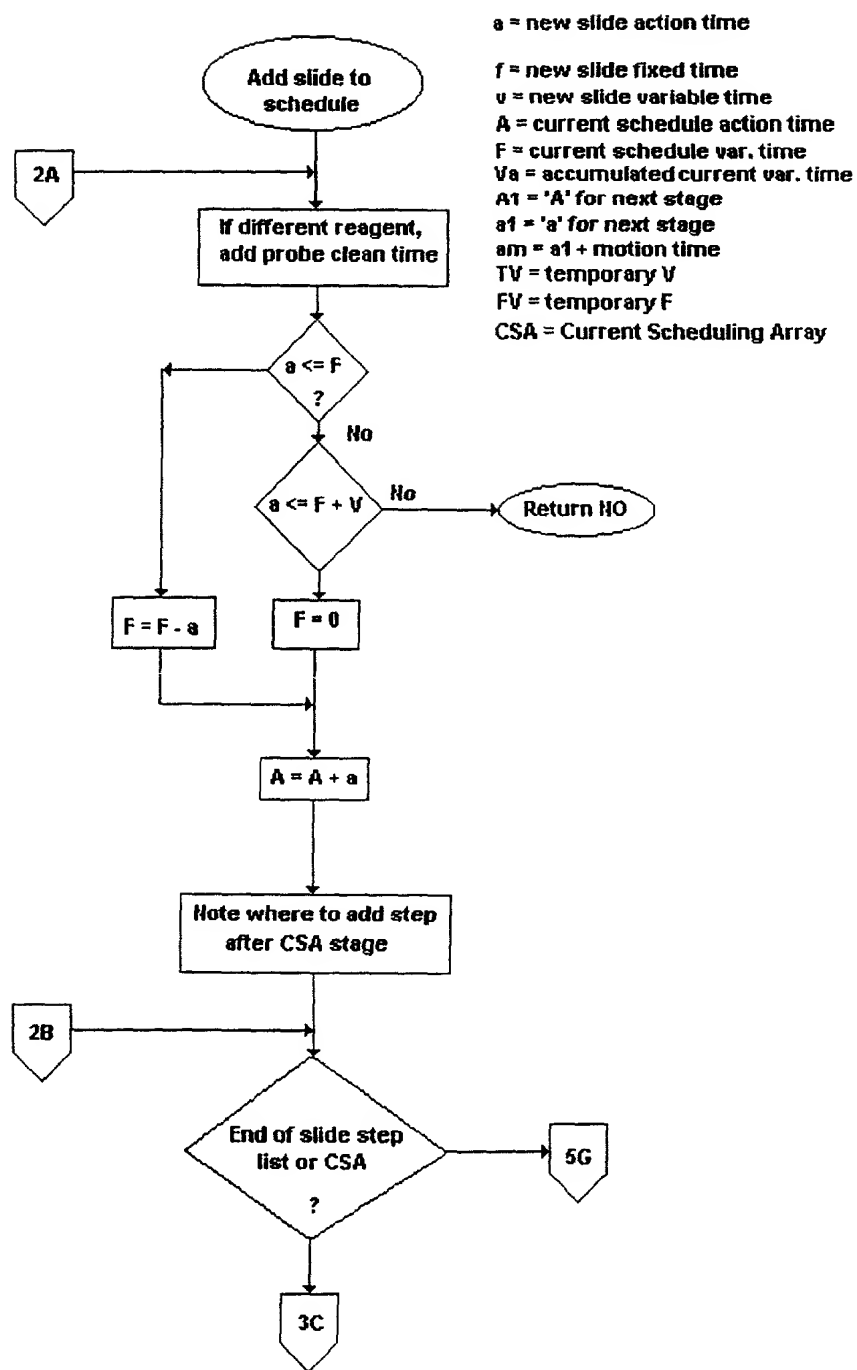


Figure 9A



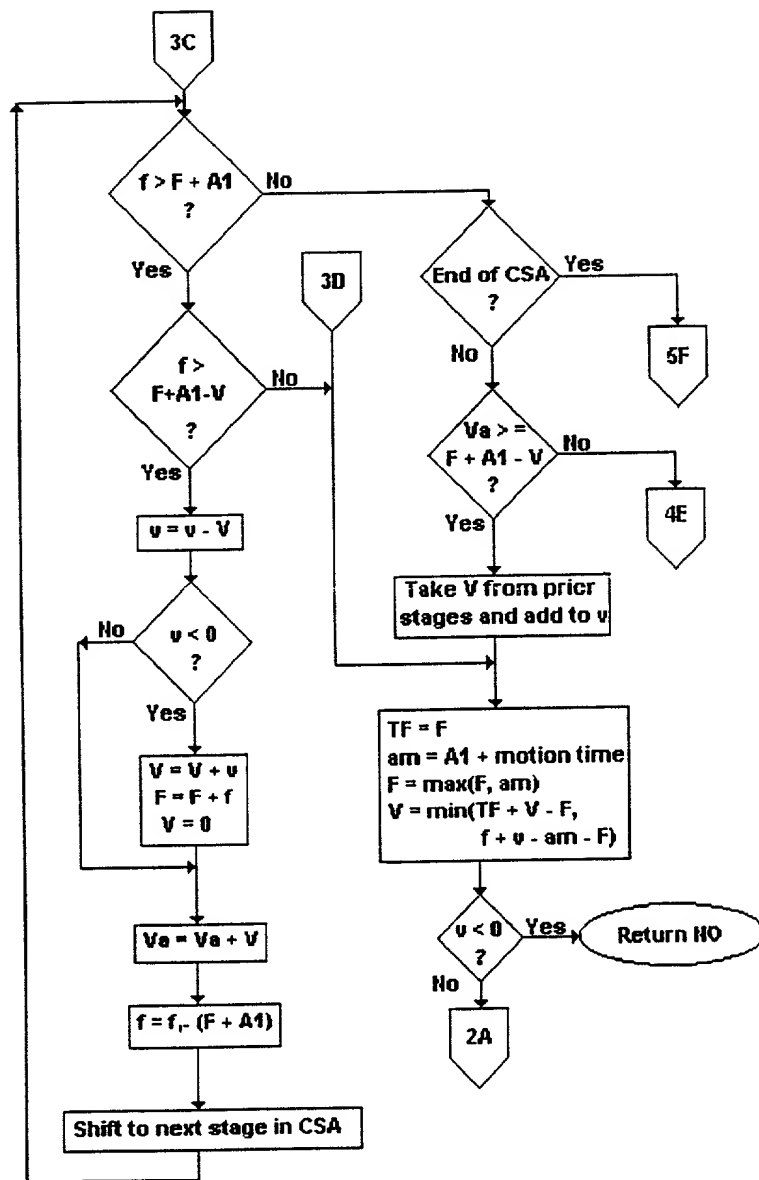


Figure 9B

Figure 9C

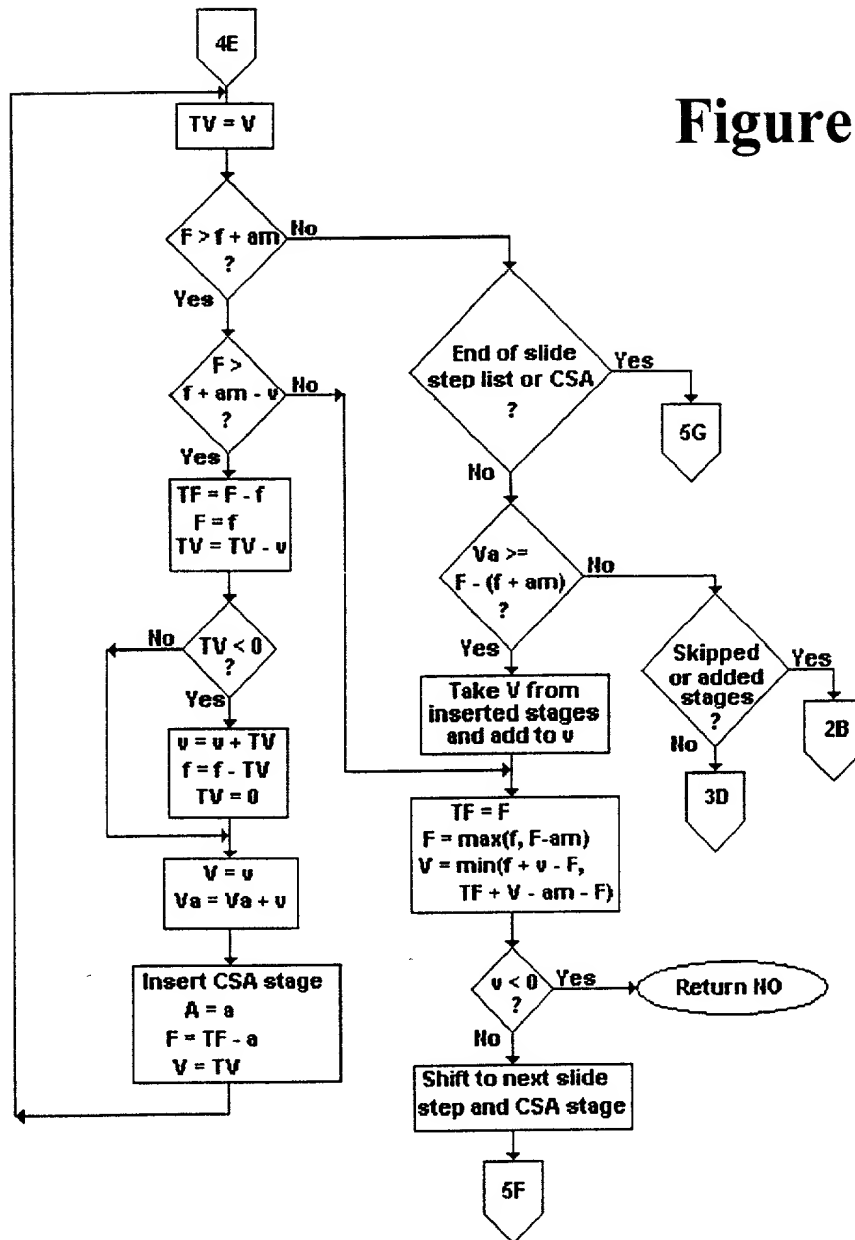

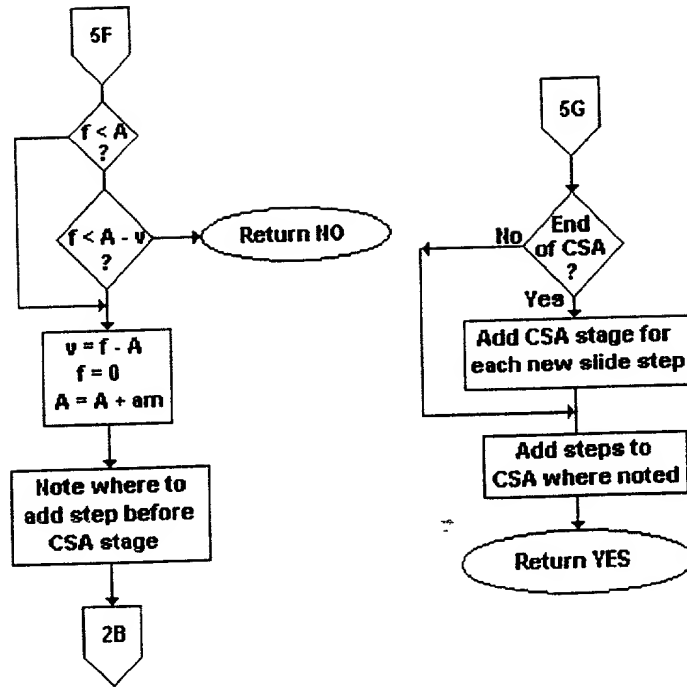

**Figure 9D**